teesa®

Massage gun

TSA0506

RELAX TOUCH **M500**

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CS

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DE

Owner's manual

EN

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FR

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Thank you for purchasing TEESA appliance. Read this instruction manual and keep it for future reference. Producer does not take responsibility for damages caused by inappropriate handling and use of the product.

SAFETY INSTRUCTIONS

In order to reduce the injuries or damages, follow basic safety precautions applied when using any electrical device, including the following:

- Read the entire instruction manual before you start installation and assembly and keep it for future reference.
- Appliance for household use only.
- Clean and store this device in accordance to instructions listed in Cleaning and storage section.
- ALWAYS turn off the device and disconnect it from the mains supply:
 - if it's not operating correctly
 - if there's an uncommon voice when using
 - before cleaning
- This appliance may be used by children who are above 8 years old and by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, if they are supervised and guided by a person who is responsible for their safety in a cautious manner for all the safety precautions being understood and followed. Children should be supervised to ensure that they do not play with the appliance. Children should not perform cleaning and servicing of the device unsupervised.
- Keep the device out of the reach of children.
- Do not use the device if it does not work correctly. Do not attempt to repair this device yourself. Only authorized and qualified personnel may repair this device. The device has no user serviceable parts inside. In case of damage, contact authorized service point for check-up or repair.
- Do not use the product for purposes different from displayed in this operating manual.
- Use only authorized accessories.
- Do not operate the appliance with wet or moist hands.
- Do not leave the appliance unattended when it's working.

- This device is not recommended to directly massage the skin to avoid possible irritation or discomfort. Massage body over clothing or place a clean and dry towel on massage part, gently press and move the appliance.
- Use the device only on soft tissues with speed setting and amount of pressure that does not produce pain or discomfort.
- Massaging the same part of body for a long time will cause bruises. Do not massage one part of your body for more than 2 minutes.
- Please keep your fingers, hair or any other body parts away from the extendable part of the device to prevent pinching.
- Before switching on, make sure massage head is not blocked.
- Prolonged continuous operation may cause the motor surface temperature to rise. Please operate carefully.


- Do NOT drop or strike the product.
- Do not use on face, head or any bony area of the body.
- Do not massage directly into the spine.
- Do not use on areas where there is bruising, pain or any discomfort.

CAUTION. Do not use this device or consult your physician if you have any concerns about your health, existing medical conditions or injuries.

4 REPLACEABLE MASSAGE HEADS

Massage head installation: Push the massage head firmly into the connection point. Make sure it has locked in place. Remove head by firmly pulling it off the gun.

CAUTION. Do not start the device without massage head in place.

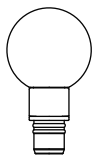
BALL HEAD

Most universal head can be used for large muscle groups such as quads, hamstrings, traps and upper back

FLAT HEAD

Use it for massaging large muscle groups such as calves and hamstrings or buttocks

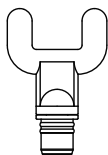FORK HEAD

Use it for massaging smaller and slimmer muscle regions such as the Achilles or along the spine

BULLET HEAD

More precise head can be used for deep tissues and joints, palms and feet

POWER ON/OFF

Press and hold the power button for 3 seconds to turn on the device to standby mode. Blue LED flashes to indicate battery level:

- Four blue lights flashing – 100%
- Three blue lights flashing – 75%
- Two blue lights flashing – 50%
- One blue light flashing – 25%

When red light starts flashing, immediately charge the device.

SELECT SPEED SETTING

When the device is on, press the power button to start the device. Blue LED is on. Keep pressing the power button to switch between four available speed levels:

- One blue light is on – first speed setting
- Two blue lights are on – second speed setting
- Three blue lights are on – third speed setting
- Four blue lights are on – fourth speed setting

After fourth press, the device will return to standby mode with LED flashing indicating battery level.

Notes:

- The machine will shut down automatically when it works continuously for 15 minutes. It is recommended to use the machine after a short rest.
- When the device is in standby mode, it will automatically shut down after 5 minutes.

CHARGING

Connect USB-C plug of the included cable to the USB-C port of the device. Connect the other plug to the AC adapter and then to the wall outlet. During charging LED light flashes.

When:

- One LED light is on – the battery is charged in 25%
- Two LED lights are on – the battery is charged in 50%
- Three LED lights are on – the battery is charged in 75%
- Four LED lights are on – the battery is charged in 100%

Please remember to disconnect the device from the power supply when the battery is fully charged (four LED lights are on).

CLEANING AND STORAGE

Clean the product with a soft dry cloth.
Store the product in a dry and cool place.

SPECIFICATION

4 speed settings
LED battery and mode indicator
15-minute auto shutoff
Compact and ergonomic design
Quiet operation

Technical data

Power: 25 W
Speed: 1600 - 3000 spm (26,6 - 50 Hz)
Operation time: 3 - 7 hours (depending on setting)
Charging time: up to 4 hours
Amplitude: up to 6 mm
Noise level: 60 dB
Material: ABS + aluminum
Battery: 2000 mAh
Power supply: 5 V / 2 A
In set: 4 massage heads, USB-C charging cable, user's manual

English
Correct Disposal of This Product
(Waste Electrical & Electronic Equipment)

(Applicable in the European Union and other European countries with separate collection systems) This marking shown on the product or its literature, indicates that it should not be disposed with other household wastes at the end of its working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources. Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can take this item for environmentally safe recycling. Business users should contact their supplier and check the terms and conditions of the purchase contract. This product should not be mixed with other commercial wastes for disposal.

Made in China for Lechpol Electronics Leszek Sp.k., ul. Garwolińska 1, 08-400 Miętno.

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