Rebel

ACTIVE


**MASSAGE
GUN**

**OWNER'S
MANUAL**

DE EN PL RO

model: RBA-6004

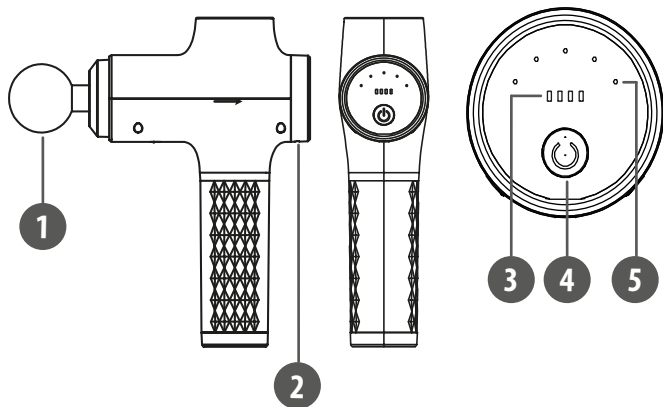6

7

8

9

10

11

12

13

	DE	EN	PL	RO
1	Massagekopf (Austauschbar)	Massage head (replaceable)	Głowica masująca (wymierna)	Cap de masaj (înlocuibil)
2	USB-C Ladebuchse	USB-C charging port	Port Țădowania USB-C	Port de Țăncăre USB-C
3	Batterieanzeige	Battery indicator	Wskaźnik poziomu baterii	Indicator baterie
4	Taste Ein/ Aus/ Geschwindigkeit	ON/OFF button/ speed switch	Przycisk ON/OFF, zmiany prędkości	Buton de alimentare/ comutator de viteză
5	Geschwindigkeitsspiegel-Anzeige	Speed level indicator	Wskaźnik poziomu prędkości	Indicator de viteză
6	Ballkopf	Ball head	Okrągła	Cap bilă
7	Gabelkopf	Fork head	Widelki	Cap furcă
8	Kugelkopf	Bullet head	Pocisk	Cap tip glonț
9	Flacher Kopf	Flat head	Płaska	Cap plat
10	Schaufelkopf	Shovel head	Łopátka	Cap tip lopátă
11	Daumenkopf	Thumb head	Kciuk	Cap tip deget mare
12	Kammkopf	Comb head	Grzebień	Cap tip pieptene
13	Bogenkopf	Arc head	Łuk	Cap tip arc

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SAFETY INSTRUCTIONS

In order to reduce the injuries or damages, follow basic safety precautions applied when using any electrical device, including the following:

- Read the entire instruction manual before you start installation and assembly and keep it for future reference.
- Appliance for household use only.
- Clean and store this device in accordance to instructions listed in Cleaning and storage section.
- ALWAYS turn off the device and disconnect it from the mains supply:
 - if it's not operating correctly
 - if there's an uncommon voice when using
 - before cleaning
- This appliance may be used by children who are above 8 years old and by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, if they are supervised and guided by a person who is responsible for their safety in a cautious manner for all the safety precautions being understood and followed. Children should be supervised to ensure that they do not play with the appliance. Children should not perform cleaning and servicing of the device unsupervised.
- Keep the device out of the reach of children.
- Do not use the device if it does not work correctly. Do not attempt to repair this device yourself. Only authorized and qualified personnel may repair this device. The device has no user serviceable parts inside. In case of damage, contact authorized service point for check-up or repair.
- Do not use the product for purposes different from displayed in this operating manual.
- Use only authorized accessories.
- Do not operate the appliance with wet or moist hands.
- Do not leave the appliance unattended when it's working.
- This device is not recommended to directly massage the skin to avoid possible irritation or discomfort. Massage body over clothing or place a clean and dry towel on massage part, gently press and move the appliance.
- Use the device only on soft tissues with speed setting and amount of pressure that does not produce pain or discomfort.
- Massaging the same part of body for a long time will cause bruises. Do not massage one part of your body for more than 2 minutes.
- Please keep your fingers, hair or any other body parts away from the extendable part of the device to prevent pinching.

- Before switching on, make sure massage head is not blocked.
- Prolonged continuous operation may cause the motor surface temperature to rise. Please operate carefully.
- DO NOT drop or strike the product.
- Do not use on face, head or any bony area of the body.
- Do not massage directly into the spine.
- Do not use on areas where there is bruising, pain or any discomfort.

CAUTION. Do not use this device or consult your physician if you have any concerns about your health, existing medical conditions or injuries.

8 REPLACEABLE MASSAGE HEADS

Massage head installation: Push the massage head firmly into the connection point. Make sure it has locked in place. Remove head by firmly pulling it off the gun.

CAUTION. Do not start the device without massage head in place.

BALL HEAD

Most universal head can be used for large muscle groups such as quads, hamstrings, traps and upper back.

FLAT HEAD

Use it for massaging large muscle groups such as calves and hamstrings or buttocks.

FORK HEAD

Use it for massaging smaller and slimmer muscle regions such as the Achilles or along the spine.

BULLET HEAD

More precise head can be used for deep tissues and joints, palms and feet.

SHOVEL HEAD

Dense muscle groups and ligaments.

COMB HEAD

Simulate finger joint.

THUMB HEAD

Palms, sole and other acupoints.

ARC HEAD

Suitable for arms, hips, thighs, calves and other large muscle groups.

POWER ON/OFF

Press and hold the power button for 3 seconds to turn on the device. Speed level and battery indicators light up. There are four LED battery indicators:

- Four LEDs are on – 100%
- Three LEDs are on – 75%
- Two LEDs are on – 50%
- One LED is on – 25% (When flashing, immediately charge the device.)

SELECT SPEED SETTING

Keep pressing the ON/OFF button to switch between four available speed levels:

Number of LEDs	Mode
One LED is on	First speed mode
Two LEDs are on	Second speed mode
Three LEDs are on	Third speed mode
Four LEDs are on	Fourth speed mode
Five LEDs are on	Fifth speed mode
Five LEDs flash alternately	Speed is set automatically according to the user's pressure

Note: The appliance will shut down automatically when it works continuously for 10 minutes. It is recommended to use the appliance after a short rest.

BLOCKING FUNCTION

When the device is operating, press the head firmly with one hand and keep it in that position for about 5 seconds to block the device - motor stops working, the device enters standby mode (LED battery and mode indicators are on). To resume operation, press the ON/OFF button. After 5 minutes in standby mode, the device automatically turns off.

TIMER

The device can operate continuously for 10 minutes, no matter if it is working or in standby mode. After the device has been blocked, it enters standby mode for maximally 5 minutes.

Example:

If the device works for 6 minutes, standby time reduces to 4 minutes what in total gives 10 minutes of operation time.

If the device works for 7 minutes, standby time reduces to 3 minutes what in total gives 10 minutes of operation time etc.

If the device operates up to 5 minutes, standby time is always 5 minutes.

Example:

If the device works for 3 minutes, standby time is 5 minutes.

If the device works for 4 minutes, standby time is 5 minutes etc.

CHARGING

Connect USB-C plug of the included cable to the USB-C port of the device. Connect the other plug to the AC adapter and then to the wall outlet. During charging LED light flashes. When:

- One LED is on – the battery is charged in 25%
- Two LEDs are on – the battery is charged in 50%
- Three LEDs are on – the battery is charged in 75%
- Four LEDs are on – the battery is charged in 100%

Please remember to disconnect the device from the power supply when the battery is fully charged (four LED lights are on).

CLEANING AND STORAGE

- Clean the product with a soft dry cloth.
- Store the product in a dry and cool place.

SPECIFICATION

6 speed settings
LED battery and mode indicator
10-minute auto shutdown
Compact and ergonomic design
Quiet operation

Technical data

Power: 60 W
Speed: 1800 - 3000 spm (30 - 50 Hz)
Operation time: up to 4 hours (depending on setting)
Charging time: up to 3,5hours
Amplitude: up to 6 mm
Noise level: 50 dB
Material: PC + ABS
Battery: 18000 mAh
Power supply: 5 V / 2 A
In set: 8 massage heads, USB-C charging cable

English
Correct Disposal of This Product
(Waste Electrical & Electronic Equipment)

(Applicable in the European Union and other European countries with separate collection systems) This marking shown on the product or its literature, indicates that it should not be disposed with other household wastes at the end of its working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources. Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can take this item for environmentally safe recycling. Business users should contact their supplier and check the terms and conditions of the purchase contract. This product should not be mixed with other commercial wastes for disposal.

Made in China for Lechpol Electronics Leszek Sp.k., ul. Garwolińska 1, 08-400 Miętno.

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