***Rebel***

**ACTIVE**

**AIR COMPRESSION LEG  
MASSAGER**

**OWNER'S  
MANUAL**

**DE**

**EN**

**PL**

**RO**

model: RBA-6007

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
## SAFETY INSTRUCTIONS

1. This device is intended for usage on legs only. Do not use it on other parts of body.
2. Before cleaning turn off the device and disconnect it from charging.
3. Clean this device after each use.
4. The device is not suitable for cleaning in a washing machine. Clean this device with soft, slightly damp cloth, without any abrasive agents.
5. If you have any doubts about using the device, consult a doctor before using it.
6. This appliance may be used by children who are above 8 years old and by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, if they are supervised and guided by a person who is responsible for their safety in a cautious manner for all the safety precautions being understood and followed. Children should be supervised to ensure that they do not play with the appliance. Children should not perform cleaning and servicing of the device unsupervised.
7. The device is not waterproof. Do not expose the device to water or moisture.
8. Store this device in original packaging, in dry place with room temperature.
9. Battery of this device is not replaceable.

## PRODUCT DESCRIPTION

1. Charging socket
2. Battery indicators
3. Mode selection button
4. Intensity selection button
5. Timer button
6. Start/stop button
7. Power button
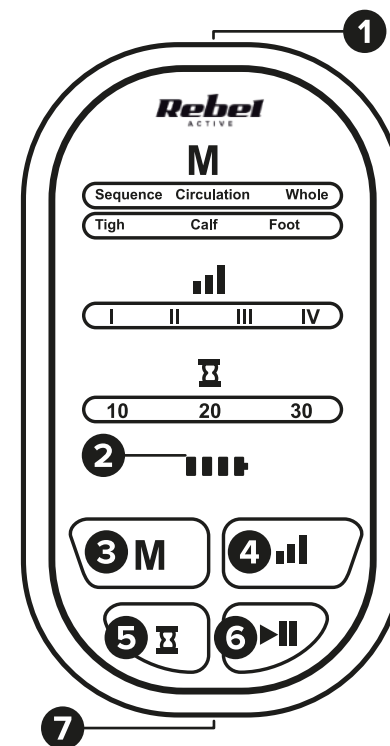
## OPERATION

- Press and hold the power button to turn the device on or off. The device will remember settings from last usage.
- Use the mode selection button to select working mode (Sequence- every part in order, circulation- circulate, whole- every part at the same time, tight, calf, foot).
- Use the intensity selection button to select one of four levels of intensity.
- Use the timer button to select the working time (in minutes).
- Use the start/stop button to start or stop the cycle.

### Charging

- Battery level is indicated by the battery level indicator. If the battery level is low, plug the charger to the device.
- Connect the charger to the USB-C charging port of the device. After finished charging, all of the indicators will glow. Disconnect the charger from the device.

## SPECIFICATION

### MAIN PARAMETERS

- Intended for: legs
- 6 working modes
- Pressure adjustment
- Timer function

### TECHNICAL DATA

- Power: 10 W
- Pressure: 50/100/150/200 mmHg
- Timer: 10/20/30 minutes
- Noise level: <55 dB
- Battery: 3000 mAh
- Charging: 5 V / 2 A
- Working time: approx. 2 hours
- Charging time: approx. 3 hours
- Charging socket: USB-C
- Dimensions: 90 x 40 cm
- Weight: 1,6 kg

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
(Applicable in the European Union and other European countries with separate collection systems) This marking shown on the product or its literature, indicates that it should not be disposed with other household wastes at the end of its working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources. Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can take this item for environmentally safe recycling. Business users should contact their supplier and check the terms and conditions of the purchase contract. This product should not be mixed with other commercial wastes for disposal.

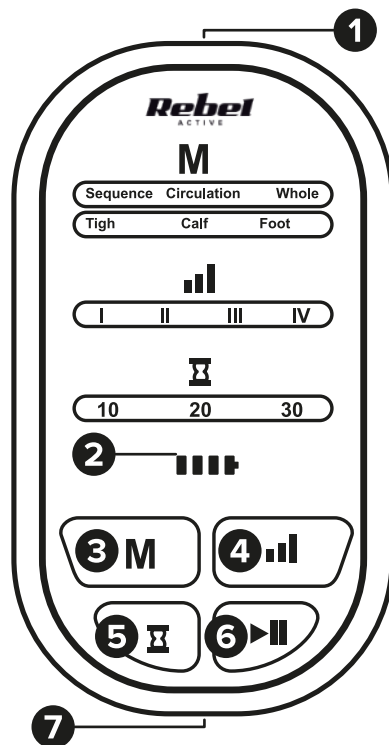Made in China for Lechpol Electronics Leszek Sp.k., ul. Garwolińska 1, 08-400 Miętno.

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