Rebel

ACTIVE

**MASSAGING BELT FOR
STOMACH**

**OWNER'S
MANUAL**

DE

EN

PL

RO

model: RBA-6008

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
SAFETY INSTRUCTIONS

1. This device is intended for use only on abdomen. Do not use it on other parts of body.
2. Before cleaning turn off the device and disconnect it from charging.
3. Clean this device after each use.
4. The device is not suitable for: pregnant women, people with irritated skin, people with heart conditions or other related illness, and people that should not be using this kind of device.
5. If you have any doubts about using the device, consult a doctor before using it.
6. This appliance may be used by children who are above 8 years old and by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, if they are supervised and guided by a person who is responsible for their safety in a cautious manner for all the safety precautions being understood and followed. Children should be supervised to ensure that they do not play with the appliance. Children should not perform cleaning and servicing of the device unsupervised.
7. The device is not waterproof. Do not expose the device to water or moisture.
8. Store this device in original packaging, in dry place with room temperature.
9. Battery of this device is not replaceable.
10. If the device is not going to be used for a longer period of time, take out the battery from the remote.

PRODUCT DESCRIPTION


Remote control

1. Massage button
2. Vibrations button
3. Intensity control
4. EMS button
5. Heating button

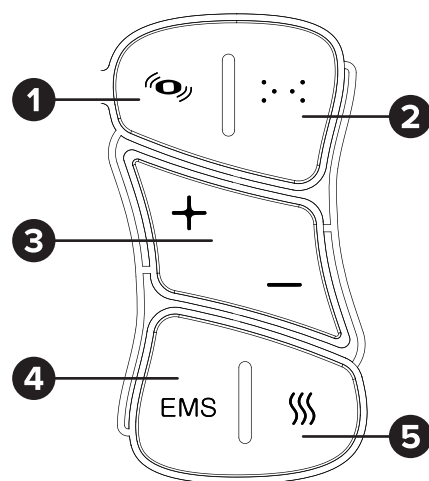Massager

1. Power button/LED indicator
2. Charging socket

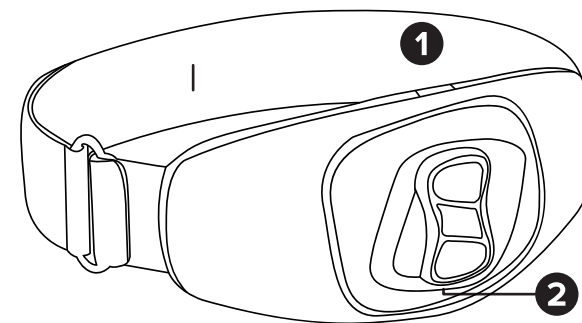OPERATION

Tuning on/off/charging

- Press and hold the power button to turn the device on or off.
- If the battery level is low, the LED indicator will glow blue. Charging will be indicated by red LED indicator. After finished charging, the LED will glow blue.

Putting on the massager

Place the massager on your abdomen area and secure it with the belt. Fasten the belt tightly enough for the massager to stay in its position, but not so tight that it results in discomfort.

Mode selection

- When the device is first powered up, it activates the first massage mode.
- To activate the massage mode, press the massage mode button. Each press will activate next massage mode, until the mode is deactivated.
- To activate the vibration mode, press the vibration mode button. Each press will activate next vibration mode, until the mode is deactivated. Massage and vibration modes can work simultaneously.
- To activate the electro stimulation mode, press the EMS button. Each press will activate next electro stimulation mode, until the

- mode is deactivated.
- To switch on or off the heating function, press the heating button. Due to increased battery consumption, it is not possible to operate massage and heating modes simultaneously.
 - The device will shut down after 15 minutes of using.

Massage intensity

To change intensity of the massage, use buttons + and -. Long beep indicates that the highest or lowest intensity level has been reached.

SPECIFICATION

MAIN PARAMETERS

- Intended for: stomach
- Massage type: vibrating, EMS
- 5 massage modes
- 4 vibration modes
- 6 EMS modes
- Heating function
- Massage intensity level adjustment
- Remote controlled

TECHNICAL DATA

- Power: 10 W
- Belt length: 95 cm
- Heating temperature: 46°C
- EMS frequency:
- Battery: 1500 mAh
- Charging: 5 V / 2 A
- Charging socket: USB-C
- Dimensions: 265 x 130 x 90 mm
- Weight: 940 g

(Applicable in the European Union and other European countries with separate collection systems) This marking shown on the product or its literature, indicates that it should not be disposed with other household wastes at the end of its working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources. Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can take this item for environmentally safe recycling. Business users should contact their supplier and check the terms and conditions of the purchase contract. This product should not be mixed with other commercial wastes for disposal.

Made in China for Lechpol Electronics Leszek Sp.k., ul. Garwolińska 1, 08-400 Miętno.

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