Rebel

ACTIVE


**NECK
MASSAGER**

**OWNER'S
MANUAL**

DE EN PL RO

model: RBA-6005

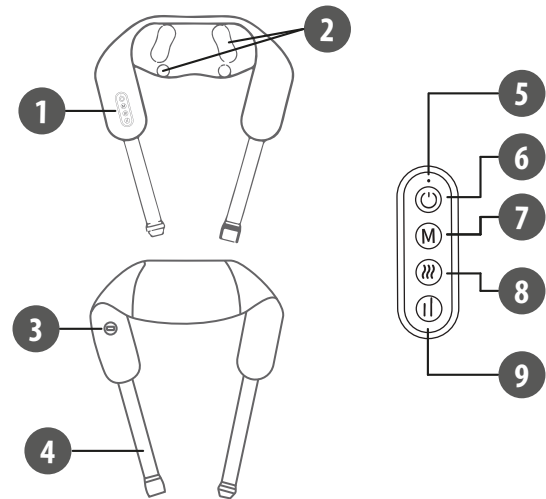	DE	EN	PL	RO
1	Kontrollpaneel	Control panel	Panel sterowania	Panou de control
2	Massageköpfe	Massage heads	Główce do masażu	Capete de masaj
3	USB-C Ladebuchse	USB-C charging port	Port ładowania USB-C	Port de încărcare USB-C
4	Rückengurt	Back strap	Zapięcie na plecy	Curea spată
5	LED Anzeige	LED indicator	Wskaźnik LED	Indicator LED
6	Taste EIN/AUS [ON/OFF]	ON/OFF button	Przycisk ON/OFF	Buton ON/OFF
7	Taste Modus	Mode button	Przycisk zmiany trybu	Buton mod
8	Taste Heizfunktion	Heating function button	Przycisk funkcji grzania	Buton funcție de încălzire
9	Taste Massageintensität	Massage intensity button	Przycisk intensywności masażu	Buton intensitate masaj

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SAFETY INSTRUCTIONS

In order to reduce the injuries or damages, follow basic safety precautions applied when using any electrical device, including the following:

- Read the entire instruction manual before you start installation and assembly and keep it for future reference.
- Appliance for household and indoor use only.
- Clean and store this device in accordance to instructions listed in „Cleaning and storage“ section.
- ALWAYS turn off the device and disconnect it from the mains supply:
 - if it's not operating correctly
 - if there's an uncommon voice when using
 - before cleaning
- This appliance may be used by children who are above 8 years old and by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, if they are supervised and guided by a person who is responsible for their safety in a cautious manner for all the safety precautions being understood and followed. Children should be supervised to ensure that they do not play with the appliance. Children should not perform cleaning and servicing of the device unsupervised.
- Keep the device out of the reach of children.
- Do not use the device if it does not work correctly. Do not attempt to repair this device yourself. Only authorized and qualified personnel may repair this device. The device has no user serviceable parts inside. In case of damage, contact authorized service point for check-up or repair.
- Do not use the product for purposes different from displayed in this operating manual.
- Use only authorized accessories.
- Do not operate the appliance with wet or moist hands.
- Do not leave the appliance unattended when it's working.
- Motor is located inside the leather cover of the massager. Do not unzip the cover and reach inside to avoid motor damage.

INTENDED USE

This device is designed for massaging the neck, shoulders, back and thighs of human body. It cannot replace medical treatment. Do not use the massager if one or several of the following applies to you. If you are unsure whether the massager is suitable for you, please consult your doctor.

Do not use the massager:

- continuously for a long time to prevent muscle pain caused by excessive massage,
- when lying flat or half lying down,
- when driving,
- if you have undergone surgery on the massaged part or have skin damage, redness, inflammation, congestion and other symptoms,
- if you suffer from heart disease, acute disease, skin disease, malignant tumor, have bleeding tendency and other blood diseases, osteoporosis, soft tissue diseases, stroke and other physical abnormalities,
- if you are pregnant woman or woman during menstrual period. If you need to use it, please follow the doctor's advice.

CHARGING

Note: Charge the device fully before first use.

1. Connect the USB-C plug of the included USB cable to the USB-C charging port on the device. Connect the other end to the AC charger and then to the power outlet.
2. During charging red LED flashes.
3. When the battery is fully charged, green LED is on. Disconnect the device from power supply.


Notes:

- When the device is charging it cannot be turned on.
- When the battery is low, the red LED indicator flashes. Please charge the device immediately.

BEFORE USE

- Before use, please remove jewelry from the massaged part or other items that may affect the massage.
- Adopt a comfortable posture and completely relax your body.
- Place the massager on the neck and shoulders* and adjust it to the most comfortable position. Fasten the strap on your back. Slip your hands into the handles.

* You can also use the massager on back and thighs.

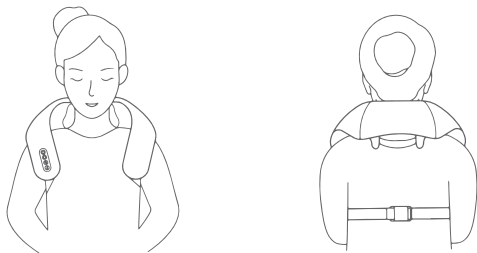OPERATION

Turn on/off: To turn on the device, press and hold the  button for 2 seconds.

By default heating function is on (red LED is on) and the intensity is set to low.

Turn off the device by pressing and holding the same button for 2 seconds. The LED goes off.

Switch mode: Press the  button to switch between two massage modes.

Heating function (on by default): Press the  button to turn on/off heating function. Red LED is on when the heating is on. When the heating function is off, green LED is on.

Set the intensity level: Press the  button to set the intensity level. There are two available settings (low and high).

Notes:

- The massage should feel pleasant and relaxing at all times. If you feel pain or discomfort on the skin, stop using it or change the contact pressure using the handles. If needed, seek medical attention immediately.
- Use as needed, within 10 minutes at a time. It is recommended to take a proper rest before another use. If used excessively, it may cause physical discomfort and skin injury. Therefore, when you use it for the first time, you need to observe your physical condition and skin changes and try to shorten the use time as much as possible.
- There is a heating surface in the massage area. Please remain cautious.

AUTO SHUT-OFF

The massager features an automatic switch-off function, which is set to the maximum operating duration of 10 minutes. If you need to continue the massage, turn on the device again.

CLEANING AND MAINTENANCE

- Before cleaning make sure the device is disconnected from power supply.
- The massage head cover can be removed and washed in the washing machine at a temperature of 30°C.
- The leather cover should be cleaned with a damp cloth. Then wipe it dry.
- Do not use any abrasives or chemical agents to clean this device.
- Store the device in a dry, cool place, out of the reach of children.

TROUBLESHOOTING

Problem	Possible cause	Possible solution
The device cannot power on	Low battery	Charge the device
The device doesn't work	Low battery	Charge the device
	The device has been working too long, overheating protection has been activated	Let the motor cool down, afterwards turn on the device
The device stops during use	Auto shut-off function has been activated	Turn on the device again
	Low battery	Charge the device
Weak massage	Low battery	Charge the device
Heating function doesn't work	The function is not turned on	Press the  button to turn on the heating function

SPECIFICATION

MAIN FEATURES

Massage different body parts
 6 silicone massage heads for shiatsu massage
 Heating function
 2 massage modes
 2 intensity levels
 Auto shut-off after 10 minutes
 Handles for massage intensity control
 Back strap
 Removable, machine washable cover
 Overheating protection

TECHNICAL DATA

Heating temperature: up to 60°C
 Operation time: up to 80 min
 Material: PU leather, polyester
 Charging port: USB-C
 Charging time: 2,5 h
 Battery: 3,7 V / 2000 mAh
 Power: 10 W
 Power supply: 5 V / 2 A
 Weight: 1,2 kg
 Dimensions: 124 x 17 x 16,3 cm
 In set: USB-C charging cable

English
Correct Disposal of This Product
(Waste Electrical & Electronic Equipment)

(Applicable in the European Union and other European countries with separate collection systems) This marking shown on the product or its literature, indicates that it should not be disposed with other household wastes at the end of its working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources. Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can take this item for environmentally safe recycling. Business users should contact their supplier and check the terms and conditions of the purchase contract. This product should not be mixed with other commercial wastes for disposal.

Made in China for Lechpol Electronics Leszek Sp.k., ul. Garwolińska 1, 08-400 Miętno.

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