Rebel

ACTIVE


**FOOT AND LEG
MASSAGER**

**OWNER'S
MANUAL**

DE EN PL RO

model: RBA-6006

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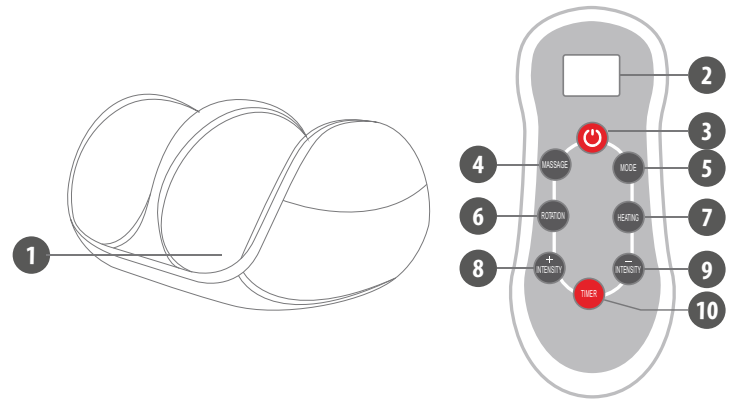SAFETY INSTRUCTIONS

In order to reduce the injuries or damages, follow basic safety precautions applied when using any electrical device, including the following:

- Read the entire instruction manual before you start installation and assembly and keep it for future reference.
- Keep the device and its cord/cable away from heat, water, moisture, sharp edges and any other factor which may damage the appliance or its cord.
- Before connecting the device to power supply socket, make sure the voltage indicated on the appliance corresponds to the voltage in the power supply socket.
- When you disconnect the power cord of the device, grasp and pull the plug, not the cord/cable.
- Device for household and indoor use only.
- Clean and store this device in accordance to instructions listed in „Cleaning and storage“ section.
- ALWAYS turn off the device and disconnect it from the mains supply:
 - if it's not operating correctly
 - if there's an uncommon voice when using
 - before disassembling the device
 - before cleaning
 - before servicing
- This appliance may be used by children who are above 8 years old and by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, if they are supervised and guided by a person who is responsible for their safety in a cautious manner for all the safety precautions being understood and followed. Children should be supervised to ensure that they do not play with the appliance. Children should not perform cleaning and servicing of the device unsupervised.
- Keep the device out of the reach of children.
- Make sure the power cord does not touch any hot or sharp elements.
- Do not use the product for purposes different from displayed in this operating manual.
- Use only authorized accessories.
- Do not operate the appliance with wet or moist hands.
- Do not leave the appliance unattended when it's working.
- Do not use the device if it does not work correctly. Do not attempt to repair this device yourself. Only authorized and qualified personnel may repair this device. The device has no user serviceable parts inside. In case of damage, contact authorized service point for check-up or repair.

- ALWAYS disconnect the product from the power mains when it is not in use.
- Motor is located inside the leather cover of the massager. Do not unzip the cover and reach inside to avoid motor damage.

INTENDED USE

This device is designed for massaging the feet, calves, thighs and forearms of human body. It cannot replace medical treatment. Do not use the massager if one or several of the following applies to you. If you are unsure whether the massager is suitable for you, please consult your doctor.

Do not use the massager:

- continuously for a long time to prevent muscle pain caused by excessive massage,
- during sleep to prevent long-term stimulation to a certain part from affecting health,
- when driving,
- if you have undergone surgery on the massaged part or have skin damage, redness, inflammation, congestion and other symptoms,
- if you suffer from heart disease, acute disease, skin disease, malignant tumor, have bleeding tendency and other blood diseases, osteoporosis, soft tissue diseases, stroke and other physical abnormalities,
- if you are pregnant or during menstrual period. If you need to use it, please follow the doctor's advice.

BEFORE USE

- Before use, please remove jewelry from the massaged part or other items that may affect the massage.
- Adopt a comfortable posture and completely relax your body.
- Place feet, calves, thighs or forearms in the massager and adjust it to the most comfortable position.

OPERATION

Connection: Connect the DC plug of the AC adapter to the DC socket of the device; connect the other end to the power outlet. OFF flashes on the display.

Turn on/off: To turn on the device, press the ☺ button. Massage starts. By default heating function is off, the intensity is set to medium, timer is set to 10 minutes.

Turn off the device by pressing the same button. OFF flashes on the display and then the display goes off.

Massage: When the device is on, press the MESSAGE button to turn on the massage (C1 appears on the display). Press the same button again to turn off the massage (CO appears on the display).

Mode: Keep pressing the MODE button to set one of three massage modes:

- AU – massage heads rotate alternately, inwards and outwards
- A1 – massage heads rotate inwards
- A2 – massage heads rotate outwards

Rotation direction (only in the massage mode): Press the ROTATION button to set the rotation direction: C1 – heads rotate inwards, C2 – heads rotate outwards.

Heating function: Keep pressing the HEATING button to set the heating level (off by default). There are three available levels- E1 (low temperature), E2 (medium temperature), E3 (high temperature). When E0 appears on the display, the heating function is off.

Intensity level: Press the INTENSITY +/- buttons to set the intensity level. There are three available settings- L1 (low intensity), L2 (medium intensity), L3 (high intensity).

Timer: Keep pressing the TIMER button to set the timer. There are three available settings: 10, 20 and 30 minutes. Set time will appear on the display. The device will operate for the time set and then it will automatically turn off.

Notes:

- The massage should feel pleasant and relaxing at all times. If you feel pain or discomfort on the skin, stop using it or change the contact pressure using the handles. If needed, seek medical attention immediately.
- Use as needed, within max. 30 minutes at a time. It is recommended to take a proper rest before another use. If used excessively, it may cause physical discomfort and skin injury. Therefore, when you use it for the first time, you need to observe your physical condition and skin changes and try to shorten the use time as much as possible.
- There is a heating surface in the massage area. Please remain cautious.

CLEANING AND MAINTENANCE

- Before cleaning make sure the device is disconnected from power supply.
- The cover is not removable. Therefore clean it with slightly damp cloth. Then wipe it dry.
- Do not use any abrasives or chemical agents to clean this device.
- When not in use, disconnect the device from power source and place it in a dry and cool place, out of the reach of children.
- During storing, do not wrap the power cord around the device.

SPECIFICATION

MAIN FEATURES

Massage different body parts
20 massage heads for shiatsu massage
Adjustable massage head rotation
3 heating settings
4 massage modes
3 intensity levels
3 timer settings

TECHNICAL DATA

Heating temperature: up to 42°C
Material: PU leather, PP
Power: 24 W
AC adapter:
 Input: 100 – 240 V; 50/60 Hz
 Output: 12 V / 2 A
Weight: 1,6 kg
Dimensions: 41 x 25,5 x 20 cm
Power cable length: 1,5 m
Remote control cable length: 1 m
In set: AC adapter

English
Correct Disposal of This Product
(Waste Electrical & Electronic Equipment)

(Applicable in the European Union and other European countries with separate collection systems) This marking shown on the product or its literature, indicates that it should not be disposed with other household wastes at the end of its working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources. Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can take this item for environmentally safe recycling. Business users should contact their supplier and check the terms and conditions of the purchase contract. This product should not be mixed with other commercial wastes for disposal.

Made in China for Lechpol Electronics Leszek Sp.k., ul. Garwolińska 1, 08-400 Miętne.

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