Rebel

ACTIVE

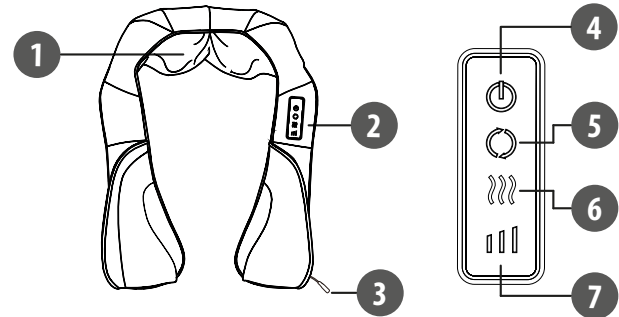
**NECK
MASSAGER**

**OWNER'S
MANUAL**

DE EN PL RO

model: RBA-6009

| | DE | EN | PL | RO |
|---|-------------------------|---------------------------|--|----------------------------|
| 1 | Massageköpfe | Massage heads | Główce do masażu | Capete de masaj |
| 2 | Kontrollpaneel | Control panel | Panel sterowania | Panou de control |
| 3 | DC Buchse | DC socket | Gniazdo DC | Mufa DC |
| 4 | Taste EIN/AUS [ON/OFF] | ON/OFF button | Przycisk ON/OFF | Buton ON/OFF |
| 5 | Taste Drehrichtung | Rotation direction button | Przycisk wyboru kierunku obrotu głowic | Buton dorecție de rotație |
| 6 | Taste Heizfunktion | Heating function button | Przycisk funkcji grzania | Buton funcție de încălzire |
| 7 | Taste Massageintensität | Massage intensity button | Przycisk intensywności masażu | Buton intensitate masaj |

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SAFETY INSTRUCTIONS

In order to reduce the injuries or damages, follow basic safety precautions applied when using any electrical device, including the following:

- Read the entire instruction manual before you start installation and assembly and keep it for future reference.
- Keep the device and its cord/cable away from heat, water, moisture, sharp edges and any other factor which may damage the appliance or its cord.
- Before connecting the device to power supply socket, make sure the voltage indicated on the appliance corresponds to the voltage in the power supply socket.
- When you disconnect the power cord of the device, grasp and pull the plug, not the cord/cable.
- Device for household and indoor use only.
- Clean and store this device in accordance to instructions listed in „Cleaning and storage“ section.
- ALWAYS turn off the device and disconnect it from the mains supply:
 - if it's not operating correctly
 - if there's an uncommon voice when using
 - before disassembling the device
 - before cleaning
 - before servicing
- This appliance may be used by children who are above 8 years old and by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, if they are supervised and guided by a person who is responsible for their safety in a cautious manner for all the safety precautions being understood and followed. Children should be supervised to ensure that they do not play with the appliance. Children should not perform cleaning and servicing of the device unsupervised.
- Keep the device out of the reach of children.
- Make sure the power cord does not touch any hot or sharp elements.
- Do not use the product for purposes different from displayed in this operating manual.
- Use only authorized accessories.
- Do not operate the appliance with wet or moist hands.
- Do not leave the appliance unattended when it's working.
- Do not use the device if it does not work correctly. Do not attempt to repair this device yourself. Only authorized and qualified personnel may repair this device. The device has no user serviceable parts inside. In case of damage, contact authorized service point for check-up or repair.

- ALWAYS disconnect the product from the power mains when it is not in use.
- Motor is located inside the leather cover of the massager. Do not unzip the cover and reach inside to avoid motor damage.

INTENDED USE

This device is designed for massaging the neck, shoulders, back, abdomen, thighs and calves of human body. It cannot replace medical treatment. Do not use the massager if one or several of the following applies to you. If you are unsure whether the massager is suitable for you, please consult your doctor.

Do not use the massager:

- continuously for a long time to prevent muscle pain caused by excessive massage,
- during sleep to prevent long-term stimulation to a certain part from affecting health,
- when driving,
- if you have undergone surgery on the massaged part or have skin damage, redness, inflammation, congestion and other symptoms,
- if you suffer from heart disease, acute disease, skin disease, malignant tumor, have bleeding tendency and other blood diseases, osteoporosis, soft tissue diseases, stroke and other physical abnormalities,
- if you are pregnant woman or woman during menstrual period. If you need to use it, please follow the doctor's advice.

BEFORE USE

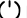
- Before use, please remove jewelry from the massaged part or other items that may affect the massage.
- Adopt a comfortable posture and completely relax your body.
- Place the massager on the neck and shoulders* and adjust it to the most comfortable position. Fasten the strap on your back. Slip your hands into the handles.

* You can also use the massager on abdomen, back, thighs and calves.

OPERATION

Connection: Connect the DC plug of the AC adapter to the DC socket of the device; connect the other end to the power outlet.

This device can be also powered from the car cigarette lighter. Use the included car charger power adapter, plug it into the DC socket on the device, then plug the car lighter plug to the cigarette lighter socket in your vehicle.

Turn on/off: To turn on the device, press the  button.

By default heating function is on (massage heads light red) and the intensity is set to medium.

Turn off the device by pressing the same button.

Rotation direction: By default the massage heads are set to automatically change direction within 1 minute. Press the  button and the massage heads will rotate inwards only. Press again and they will rotate outwards only.

Set the intensity level: Press the  button to set the intensity level. There are three available settings (low, middle and high).

Heating function (on by default): Press the  button to turn on/off heating function. Massage heads light red when the heating is on.

Notes:

- The massage should feel pleasant and relaxing at all times. If you feel pain or discomfort on the skin, stop using it or change the contact pressure using the handles. If needed, seek medical attention immediately.
- Use as needed, within 15 minutes at a time. It is recommended to take a proper rest before another use. If used excessively, it may cause physical discomfort and skin injury. Therefore, when you use it for the first time, you need to observe your physical condition and skin changes and try to shorten the use time as much as possible.
- There is a heating surface in the massage area. Please remain cautious.

AUTO SHUT-OFF

The massager features an automatic switch-off function, which is set to the maximum operating duration of 15 minutes. If you need to continue the massage, turn on the device again.

CLEANING AND MAINTENANCE

- Before cleaning make sure the device is disconnected from power supply.
- The cover is not removable. Therefore clean it with slightly damp cloth. Then wipe it dry.
- Do not use any abrasives or chemical agents to clean this device.
- Store the device in a dry, cool place, out of the reach of children.

SPECIFICATION

MAIN FEATURES

Massage different body parts
8 massage heads for shiatsu massage
Heating function
Adjustable massage head rotation
3 intensity levels
Auto shut-off after 15 minutes
Handles for massage intensity control
Overheating protection

TECHNICAL DATA

Heating temperature: up to 42°C
Material: PU, nylon
Power: 24 W
AC adapter:
Input: 100 – 240 V; 50/60 Hz
Output: 12 V / 2 A
Weight: 1,3 kg
In set: AC adapter, car charger power adapter

English
Correct Disposal of This Product
(Waste Electrical & Electronic Equipment)

(Applicable in the European Union and other European countries with separate collection systems) This marking shown on the product or its literature, indicates that it should not be disposed with other household wastes at the end of its working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources. Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can take this item for environmentally safe recycling. Business users should contact their supplier and check the terms and conditions of the purchase contract. This product should not be mixed with other commercial wastes for disposal.

Made in China for Lechpol Electronics Leszek Sp.k., ul. Garwolińska 1, 08-400 Miętno.

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