reesa®

Sonic toothbrush

TSA8082-B
TSA8082-W

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CS

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DE

Owner's manual

EN

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FR

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Thank you for purchasing TEESA appliance. Read the entire instruction manual before you start installation and assembly and keep it for future reference. Producer does not take responsibility for damages caused by inappropriate handling and use of the product.

SAFETY INSTRUCTIONS

Follow basic safety rules when operating the device, and follow instructions below:

1. Read this instruction manual carefully, even if you are familiar with the appliance.
2. Before connecting the device to power supply socket, make sure the voltage indicated on the appliance corresponds to the voltage in the power supply socket.
3. ALWAYS turn off the device and disconnect it from the mains supply:
 - if it's not operating correctly
 - if there's an uncommon voice when using
 - before disassembling the device
 - before cleaning
 - before servicing
4. When you disconnect the power cord, grasp and pull the plug, not the cord/cable.
5. This appliance may be used by children who are above 8 years old and by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, if they are supervised and guided by a person who is responsible for their safety in a cautious manner, and all the safety precautions are understood and followed. Children should not play with this device. Children should not perform cleaning and servicing of the device unsupervised.
6. Do not use this appliance if it's not operating properly or has been dropped into water. Do not attempt to repair this appliance yourself, only authorized and qualified personnel may repair this device. This product has no user serviceable parts. In case of damage, contact an authorized service point for check-up/repair.
7. Keep the device and its cord/cable away from heat, water, moisture, sharp edges and any other factor which may damage the appliance or its cord.
8. Do not use the product for purposes/uses different from what displayed in this operating manual.
9. Use only original/authorizes accessories.
10. DO NOT USE while bathing!
11. Do not place or store product where it can fall or be pulled onto a tub or sink.
12. Do not reach for a product that has fallen into water. Unplug from power supply socket immediately.
13. Do not handle plug or the appliance with wet hands.
14. Product for household use only. Not for commercial use.

IPX7

This product is water-resistant based on the **IPX7** rating, which tests submersion up to 1 meter for up to 30 minutes. Despite this classification, your device is not impervious to water damage in any situation.

The device is electrically safe and is designed for use in the bathroom- it is suitable for cleaning under running water but it is not meant to be used underneath the water.

Make sure the plug at the bottom of the toothbrush handle is always tightly sealed to prevent water from getting into the charging socket.

Intended use



This product is designed solely for cleaning teeth, gums and tongues. Do not use it for other purposes.

Initial use of this product may cause tooth discomfort or minor gum bleeding due to the gums not being accustomed to the high-efficiency cleaning. This phenomenon should cease within 1-2 weeks. If it persists beyond 2 weeks, discontinue use and consult a dentist.

Users with dental restorations, braces, or implants should use this product under the guidance of a dentist.

CHARGING

1. Connect the included USB cable to the USB-C socket located at the bottom of the toothbrush handle.
2. Connect the other plug of the USB cable to the AC charger and then to the power outlet.

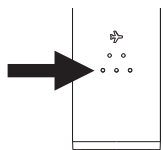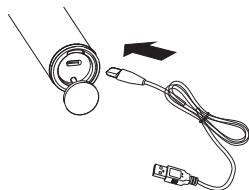3. During charging LED indicators in the lowest row flash.

4. Once fully charged, all three LEDs are on.
5. Disconnect the toothbrush from power supply.


Notes:

- Make sure the battery is fully charged before first use.
- During charging the bottom of the handle heats up. It is a normal phenomenon.
- If the battery level drops below 15%, the LEDs indicating the currently selected intensity level will flash. If the battery is close to running out, the 3 intensity level indicators will light up and the device will turn off after a while.


CONFIGURATION

3 SETTINGS

Select desired setting by pressing the **settings button**.

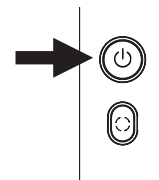
Settings indicators (blue) are those in the **middle**.

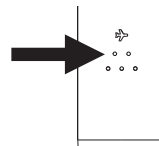1 VIBRATION (adjust the intensity of vibrations)

Indicator: The indicator on the **left** is on.


Vibration levels: 7°, 8°, 9°, 10°, 11°

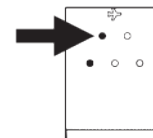2 SWEEP ANGLE (adjust the sweep angle (the angle at which the toothbrush head moves) to cover more surface area)

Indicator: The indicator on the **right** is on.


Sweep angles: 6°, 18°, 30°, 42°, 54°

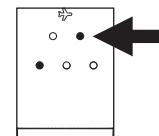3 SWEEP SPEED (regulates the speed at which the toothbrush head sweeps back and forth)


Indicator: **Both** indicators are on.

Sweep speeds: 90, 150, 210, 270, 330 times per minute

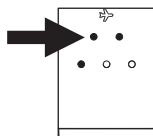
5 LEVELS

Select the desired level by pressing the **level selection button**.

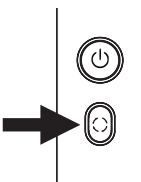Level selection indicators (white) are those at the **bottom**.

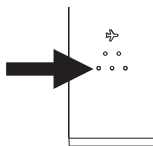First level	Second level	Third level	Fourth level	Fifth level

PERSONALIZATION

The toothbrush operation consists of: **vibration, sweep angle** and **sweep speed**. The combination of those three settings allows to achieve an individual program tailored to the user's needs.



1. Press the settings button to enter first setting - **VIBRATION**. Then keep pressing the level selection button to set one of five available vibration levels.
2. Press the settings button again to enter second setting - **SWEEP ANGLE**. Then keep pressing the level selection button to set one of five available sweep angle levels.
3. Press the settings button again to enter third setting - **SWEEP SPEED**. Then keep pressing the level selection button to set one of five available sweep speed levels.

NOTES:

- When the toothbrush is switched on, all three settings turn on simultaneously on the first level. If the user stops the configuration on one of the settings (vibration, sweep angle or sweep speed) then the two remaining settings will work on the first level or on the previously selected level if the toothbrush has been already used.
- Please note, you will only be able to change the setting within 5 seconds from the last time you pressed the button. Once the 5 seconds are up and you press the settings button, the appliance will turn off.

BRUSHING


1. Attach the toothbrush head to the shaft as shown in the diagram, leaving approximately 0,6 mm of gap between the toothbrush head and the end of the handle.
2. Rinse the toothbrush head and apply an appropriate amount of toothpaste. Position the bristles at a small angle against the gum line and place them on the teeth.
3. Apply gentle pressure and move the toothbrush head slowly and with small movements between the teeth. Avoid excessive force.
4. To clean the entire mouth more effectively, divide the mouth into four quadrants. Start with the first quadrant (outer upper teeth), brushing each quadrant for 30 seconds. Every 30 seconds, the brushing action will pause to remind you to move to the next quadrant. The two-minute timer will automatically turn off when the brushing time is complete.
5. Clean the toothbrush head and handle by rinsing them with water. Gently tap the toothbrush head to ensure no food residue remains in the bristles.
6. Remove the brush head in the direction indicated in the diagram and shake off any excess water.

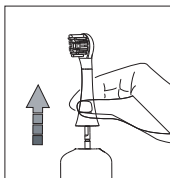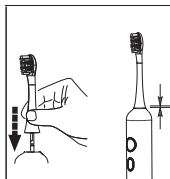
ANTI-TOUCH PROTECTION

This appliance features *anti-touch protection*. When enabled, the appliance will not react to any button press. It prevents accidental switch on for example when your toothbrush is in the suitcase.

To turn on the function press and hold the **settings button** for approx. 5 seconds.

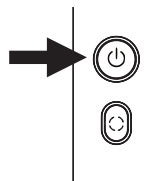The **anti-touch protection LED indicator** lights up.

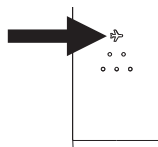To switch off, press and hold the same button for approx. 5 seconds. When the indicator goes off, it means that the function has been disabled.

LAST SETTING MEMORY

The appliance is equipped with last setting memory. It will boot in the configuration it has been used lately.

CLEANING AND MAINTENANCE

Head cleaning:

1. Gently remove the toothbrush head from the handle.
2. Rinse the toothbrush head and bristles with clean water to remove food residue and toothpaste foam.
3. It is recommended to remove the toothbrush head from the handle and rinse the connection area at least once a week to prevent bacterial growth.

Handle cleaning:

1. Wipe the handle surface with a damp cloth to remove dirt and debris.

Brush head replacement

It is recommended to replace the brush head with a new one **every three months** when using twice a day.

Replacement brush heads (soft and medium) are available for purchase online at www.teesa.pl or www.rebelelectro.com.

Toothbrush travel case is available for purchase online at www.teesa.pl or www.rebelelectro.com.

SPECIFICATION

MAIN FEATURES

Vibration level, sweep angle and sweep speed adjustment for personalized program
 31 000 strokes/min. (in each setting)
 2 min. auto-off, 30 sec. timer
 Anti-touch protection
 Last setting memory
 Charging indicator

TECHNICAL DATA

IP rating: IPX7
 Noise level: ≤65 dB
 Usage on one charge: up to 30 days
 Battery capacity: 1100 mAh
 Power supply: 5 V / 1 A
 In set: 2x soft brush head, USB cable

English
Correct Disposal of This Product
(Waste Electrical & Electronic Equipment)

(Applicable in the European Union and other European countries with separate collection systems) This marking shown on the product or its literature, indicates that it should not be disposed with other household wastes at the end of its working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources. Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can take this item for environmentally safe recycling. Business users should contact their supplier and check the terms and conditions of the purchase contract. This product should not be mixed with other commercial wastes for disposal.

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